

6970EFM Freezer/Mill®

Grinding and
Pulverizing



Tissue
Homogenization
Cell Lysis




EXTERNAL LOADING



ENCLOSED SYSTEM WITH LIQUID NITROGEN AUTOFILL



HIGH THROUGHPUT

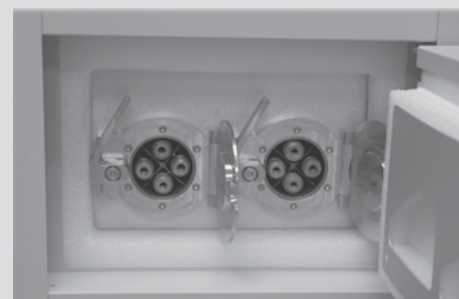

MULTIPLE VIAL SIZES



SPEX SamplePrep®

6970EFM

Freezer/Mill®

THE LATEST ADVANCES IN CRYOGENIC GRINDING

- ■ For samples that are impossible to grind at ambient temperature
- ■ Backed by SPEX SamplePrep's 40 years of experience in cryogenic grinding
- ■ Ideal for grinding bone, teeth, hair, tough plant and animal tissue, textiles, plastics, pharmaceuticals and more
- ■ DNA/RNA extraction without shearing
- ■ No loss of volatile compounds

Accuracy

- ■ Samples are maintained at Liquid Nitrogen (LN) temperatures
- ■ Vials can be cleaned or discarded to eliminate possibility of contamination

Safety and Convenience

- ■ Totally enclosed LN auto-fill system with no LN exposure
- ■ CE compliant

Durability

- ■ Unique maintenance free design, only one moving part

Versatility

- ■ Entire system is controlled by a touch screen interface capable of storing up to 10 grinding programs for rapid process replication and consistency
- ■ Grinding vials available in a range of sizes from 0.6 mL to 50 mL

FIELD OF APPLICATION

- ▶ Biotechnology
- ▶ Plant Research and Crop Science
- ▶ Pesticide Residues
- ▶ Medicine and Pharmaceuticals
- ▶ Marine Biology
- ▶ Forensics
- ▶ Polymers
- ▶ RoHS WEEE
- ▶ Food Science
- ▶ Materials Science
- ▶ Nanotechnology and Mechanical Alloying
- ▶ Glass and Ceramics
- ▶ Geology and Mining
- ▶ Environmental Science
- ▶ Cement, Slag and Fly Ash
- ▶ Catalysts
- ▶ National Defense and Homeland Security
- ▶ Veterinary Drugs



6970D

includes two grinding chambers, increasing grinding capacity to eight 6751 grinding vials. Chambers can be programmed for independent operation and different protocols, providing flexibility and speed.



6970M

includes one grinding chamber and one pre-cool/storage chamber. Can be upgraded to a 6970D.

For more information please contact 732-623-0465 or visit www.spexsampleprep.com/freezermill

ACCESSORIES

STANDARD GRINDING VIALS

- Sample integrity is maintained in tightly sealed vials
- Polycarbonate center cylinders with steel end plugs and impactors
- Polycarbonate center cylinders are reusable or can be used as disposables
- Three sizes: small samples up to 4mL, mid-size samples for up to 25mL, and large samples for up to 50mL
- Up to four small vials per grinding chamber using multi-vial adapter



6751 grinding vials

SPECIAL GRINDING VIALS

- Alternative materials available for most vial types
- Poly-Vial sets for grinding soft tissues without metal contamination
- Cr-Free Vial set for RoHS/WEEE testing of electronic components
- Micro-Vial set holds three micro-vials for samples from 0.1-0.5mL
- Steel Center Cylinders for when contact with plastic is not advised



Cr-Free Impactor

Poly Vial

VIAL EXTRACTORS AND VIAL RACKS

- Vial extractors are used to remove end plugs from center cylinder
- Free vial extractor and rack accessory package supplied with each 6970 Freezer/Mill



AUTO-EXTRACTOR™

- The 6814 Auto-Extractor is an electrically powered version of the manual vial extractors
- Increases throughput by speeding up end plug removal
- Works with large and small vials



SPEX SAMPLEPREP OFFERS THE COMPLETE SOLUTION FOR AUTOMATED CRYOGENIC GRINDING

- Portable cryogenic 50L Dewar designed for storing and dispensing liquid nitrogen through the auto-fill system
- Cryogenic Transfer Hoses to connect Dewar to the auto-fill outlet in the Freezer/Mill
- Cryogenic gloves—designed to protect the hands and arms from ultra-cold materials



BEFORE & AFTER SAMPLES



APPLICATION NOTES

Examples

- SP005 – RNA Extraction from Cartilage Tissue using Cryogenic Grinding
- SP006 - Preparation of Rice Husks for Content Analysis
- SP007 - Grinding Polymers for Qualitative and Quantitative Analyses
- SP015 - Molecular Effects of Pollution in Cetacean Samples (RNA Extraction from Whale Skin and Blubber)

Please visit

www.spexsampleprep.com/appnotes

for more FREEZER/MILL application notes.

6970EFM

Freezer/Mill®

ABSTRACTS

Biotechnology & Forensics

RNA Extraction from Knee Cartilage Tissue

A study to determine the connection between arthritic diseases and certain genes was conducted. Knee cartilage tissue was taken from patients and snap frozen in liquid nitrogen for storage. Samples were ground in the large Freezer/Mill using small vials and a multi-vial adapter. Grinding was completed in three cycles of two minutes each. The total processing time of 12 minutes included precool, grinding, and cooling between grinding cycles. RNA was subsequently purified using a phenol-chloroform extraction method with yields of RNA between 2 & 20 µg.

Polymers & Materials Research

Grinding Polymers for Qualitative and Quantitative Analysis

Compact synthetics (pellets, webbing fabric, textiles) and polymer yarns, fleece material and filaments were successfully ground in the small Freezer/Mill. Cryogenic grinding was required for accurate qualitative and quantitative analysis of additives such as antioxidants using HPLC and INOC-Test. The Freezer/Mill efficiently reduced these samples to particle sizes of less than 500 µm in less than 15 minutes without thermal degradation.

Crop Research

Preparation of Rice Husks for Content Analysis

A study to analyze the quantity of vitamins and important minerals in rice husks was performed. 8g of sample was ground in the large Freezer/Mill. Rice husks are extremely tough and in 3 grinding cycles of 3 minutes each, very good homogenization of the rice husks was achieved. The analytical fineness of the homogenized sample enabled efficient extraction of the desired components for analysis.

SPEX SamplePrep. 

Your Science is Our Passion.®

Customer Service

Sampleprep@spexsp.com

15 Liberty Street • Metuchen, NJ 08840

732-623-0465 • Fax 732-906-2492

www.spexsp.com

2010

THE ULTIMATE PLANT & ANIMAL TISSUE PULVERIZER & CELL LYSER

The 2010 Geno/Grinder

- ▶ Specifically designed for vigorous up-and-down shaking of deep-well titer plates and sample vials from 0.6 to 50mL
- ▶ Used to prepare tissue for extractions of nucleic acid, protein, and other constituents such as pesticides
- ▶ Sample material includes seeds, stems, roots, leaves and soft animal tissue
- ▶ Samples are pulverized with a large choice of grinding media, with or without buffer
- ▶ Unique vertical motion with variable rate settings for mixing and pulverizing samples

For more information on the Geno/Grinder, please visit:
www.spexsampleprep.com/genogrinder

Geno/Grinder®



SPEX is proudly represented in Australia by
AXT Pty. Ltd.
1/3 Vuko Pl., Warriewood
NSW 2102 Australia
T. +61 (0)2 9450 1359 F. +61 (0)2 9450 1365
W. www.axt.com.au E. info@axt.com.au

SPEX SamplePrep®

SPEX SamplePrep
15 Liberty Street
Metuchen, NJ 08840-9978
Tel: 732.623.0465
Fax: 732.906.2492
800.522.7739 (800.LAB.SPEX)
E-mail: sampleprep@spexcsp.com
www.spexsampleprep.com

European Division
SPEX CertiPrep Ltd
2 Dalston Gardens
Stanmore, HA7 1BQ, UK
Tel: +44 (0)20 8204 6656
Fax: +44 (0)20 8204 6654
E-mail: sales@spexcertiprep.co.uk
Web: www.spexcsp.com

DEMO PROGRAM: TRY IT BEFORE YOU BUY IT

SPEX SamplePrep gives you the opportunity to try out our sample preparation equipment in your laboratory for one week.

Best of all, our Demo Program is COMPLIMENTARY!!!

For more information please contact 1-800-LAB-SPEX or visit www.spexsampleprep.com/freezermill

Once you try our quality products you'll ask yourself how you ever lived without them.